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**HoMedics
Service Center Dept. 168
3000 Pontiac Trail
Commerce Township, MI
48390**

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Tri-Fitness™ Plus HealthStation™

Know exactly what you're made of.™

How are Body Fat (BF), Total Body Water (TBW), and Muscle Mass (MM) % Estimated?

This HealthStation uses BIA (Bioelectrical Impedance Analysis) to measure BF%, TBW% and MM%. BIA sends an unnoticeable, low intensity electrical signal through the body. The signal is very low and causes no bodily harm. The signal flows with a different degree of difficulty depending upon the individual's body fat level. This difficulty is called electrical impedance and is how BF%, TBW% and MM% can be estimated.

Why is Measuring Your Body Fat % Important?

The goal in proper weight loss is to increase lean muscle mass (which is the shape and contour of your body) and to decrease body fat for a higher burning metabolism.

There are two basic types of body fat: essential fat and storage fat.

Essential Fat is needed for normal, everyday physiological functions within the human body. It is located in and around muscle tissues, nerve cells, bone marrow, intestines, liver and lungs. Without essential fats, normal body functions begin to degrade and ultimately shut down.

Storage Fat is found and stored just below the skin's surface and around major organs. Storage fat acts as an insulator for the body, provides energy for metabolism and is a cushion against trauma. Men tend to have storage fat around the waistline and women tend to carry storage fat in the hips and upper thighs. Though some storage fat is necessary for normal body functions, excessive amounts can be detrimental to your health and overall quality of life.

Traditionally, total body weight determined whether or not a person was obese. However, weight changes can be affected by both muscle and fat gain/loss. A change in weight alone does not tell you whether it was a gain/loss of muscle or fat. In weight management, you want to maintain muscle mass while body fat is lost. Monitoring BF% is an important step toward successful weight management and body health.

This instrument is a personal monitor and should be used in a consistent manner for the most accurate reading. While readings might vary between this analyzer and other methods of measurement (hydro-densitometry or hand held calipers), your **changes** in body fat percentage shown by this instrument will be reflected accurately.

The optimal Body Fat % of an individual varies according to age and gender. Use the table below as a guide:

BODY FAT % RANGES FOR MEN[†]

Rating	Age				
	20-29	30-39	40-49	50-59	60+
Low	<13	<14	<16	<17	<18
Optimal	14-20	15-21	17-23	18-24	19-25
Moderate	21-23	22-24	24-26	25-27	26-28
High	>23	>24	>26	>27	>28

BODY FAT % RANGES FOR WOMEN[†]

Rating	Age				
	20-29	30-39	40-49	50-59	60+
Low	<19	<20	<21	<22	<23
Optimal	20-28	21-29	22-30	23-31	24-32
Moderate	29-31	30-32	31-33	32-33	33-35
High	>31	>32	>33	>34	>35

[†] Source: University of Illinois Department of Food Science and Human Nutrition. Body Fat Percentage Calculator. www.ag.uiuc.edu/~food-lab/ai/bfc.html

Why is Measuring Your Total Body Water Percentage Important?

Monitoring your level of body water can be a useful tool for health maintenance. The more lean muscle you gain as you exercise and get in shape, the lower your body fat and the higher your body water will be.

Water is an essential body component and is one of the body's overall health indicators. It makes up approximately 50-70% of the body's weight and is present proportionally more in lean tissue as compared to fat tissue.

For example:

73% of lean muscle tissue is water

80% of blood is water

25% of fat is water

Additionally, water is a medium for biochemical reactions that regulate body functions

- Removes waste products from cells for excretion in urine and sweat
- Provides form to cells which helps maintain body temperature
- Provides moisture to skin, mucous membranes and cushions vital organs

- Lubricates joints and is a component of many body fluids
- Helps flush fat (triglycerides) out of your system
- Helps you feel full and reduces your appetite
- Aids in proper digestion and helps prevent constipation
- Helps reduce the chances of developing kidney stones

The normal TBW% of an individual varies according to BF% and gender. The table below may be used as a guide:

BODY WATER CHART^{††}

	BF % Range	Optimal TBW % Range
Men	4 to 14%	70 to 63%
	15 to 21%	63 to 57%
	22 to 24%	57 to 55%
	25 and over	55 to 37%
Women	4 to 20%	70 to 58%
	21 to 29%	58 to 52%
	30 to 32%	52 to 49%
	33 and over	49 to 37%

^{††} Source: Derived from Wang & Deurenberg: "Hydration of fat-free body mass". American Journal Clin Nutr 1999; 69 833-841.

Your TBW% may vary according to your hydration level. For better accuracy, avoid changes in hydration level before taking a measurement. The accuracy in estimating TBW% will also decrease for individuals suffering from illnesses where water tends to accumulate in the body. Please remember that the scale represents only a good approximation of your total body water.

Why Should I Know My Muscle Mass %?

According to the American College of Sports Medicine (ACSM), lean muscle mass may decrease by nearly 50 percent between the ages of 20 and 90. If you don't do anything to replace that loss, you're losing muscle and increasing fat. It is also important to know your muscle mass % during weight reduction. At rest, the body burns approximately 50 additional calories for each pound of muscle gained. Some advantages of gaining muscle mass include:

- Reversing the decline in strength, bone density and muscle mass that accompanies age
- Maintenance of flexible joints
- Guide weight reduction when combined with a healthy diet

The normal muscle mass for men is greater than 40%. For women, normal muscle mass is greater than 34%.

What is BMR?

Your BMR (Basal Metabolic Rate) is the energy (measured in calories) expended by the body, at rest, to maintain normal body functions. This continual work makes up about 60-70% of the calories the body uses, and includes the beating of your heart, respiration, and the maintenance of body temperature. Your BMR is influenced by a number of factors, including age, weight, height, gender, dieting and exercise habits.

Why Should I Know My BMR?

If you know how many calories your body needs to function, you can calculate exercise and calorie consumption into your diet and weight loss program to ensure that you are eating and exercising at your proper level.

You can boost your BMR by building muscle and reducing body fat. This increases the body's ability to burn more calories while the body is at rest. When you raise your BMR, your body is less likely to store fat. The amount of increase in BMR varies from person to person, but even a modest increase in activity should cancel out the body's tendency to decrease BMR when you cut calories in your diet.

Note: Ideal body fat content, total body water, muscle mass, BMR and exercise guidelines are not the same for all people. Age, sex, heredity and overall health are variables in this measurement. Always consult your physician to determine what is most ideal for you.

Fitness Levels

This HealthStation has three fitness levels for more precise readings.

Fitness Level 1

- Normal mode
- Use at the beginning of an exercise/diet program
- Person with relatively low level of activity
- Physical activity: less than 2 times per week, less than 20 minutes per session

Fitness Level 2

- Use as you progress through an exercise/diet program when benefits of overall fitness are achieved.
- Physical activity: moderate daily exercise; 5 times per week, at least 20 minutes per session

Fitness Level 3

- Highly active individual with high endurance and strength
- Physical activity: 5 times per week, at least 60 minutes per session

When Should I Use the Body Fat, Body Water and Muscle Mass Functions?

For maximum accuracy and repeatability, it is recommended that the body fat, body water and muscle mass functions be used at approximately the *same time of day*, under the *same circumstances*, e.g. before breakfast in the morning. It is also a good practice to avoid swings in the body's hydration level prior to taking the measurement. By establishing the pattern of your body's fluctuation over a period of one month, you can determine your average level of body fat, body water and muscle mass %. Women will have more fluctuation in water content than men and need to monitor their reading over a longer time period to establish a personal average. You must determine this as a reference point before measuring progress in a health program.

You must wait several hours before taking a reading when: **drinking coffee or alcohol, taking diuretic medications or exercising. These all affect your level of hydration and the accuracy of this analyzer.**

The **condition of the skin on the bottom of your feet** can also affect the reading. The natural effects of aging or activity can make this skin hard. Take the reading with clean, slightly moist feet for best accuracy.

Caution

BIA method determines your body fat, body water and muscle mass % by sending a harmless signal through the body. Do not use this product if you have a **pacemaker, artificial limbs, contraceptive devices, metal plates/screws or other internal medical devices. When in doubt, contact your physician.**

This HealthStation will give accurate body fat, body water, and muscle mass % plus BMR for a majority of people, but it is not intended for use by the following groups:

Children: Anyone under 10 years of age
Pregnant Women

Battery Information

Your unit has a lithium battery installed that will provide a long life of operation when used in a typical household. Under heavy usage, the battery may eventually need to be replaced. Remove the battery if the product is not used for a long time.

Prior to Use

Remove red tab on bottom of unit to activate battery.
This can be done without opening battery compartment.

Getting Started

Select pounds (Lb/in), kilograms (Kg/cm) or stones (St/in) using the slide switch on the bottom of the unit.

Weight Measurement Only

You can use this function anytime-with shoes or without- by stepping on the HealthStation.

Place the analyzer on a firm, flat surface to achieve best accuracy.

Step onto the analyzer and stand still while your weight is calculated.
Your weight will be displayed on the LCD screen as shown below

The unit will turn off by itself after use.

Programming Personal Data into Memory

NOTE: In Personal Data mode, the unit returns to standby mode if no button is pressed within 30 seconds.

Press '**Set**' button; the number P-01 will appear and flash.
Use the (▲) or (▼) button to select a user number between P-01 and P-10. Press '**Set**' to confirm.


Set Height – Use the (▲) or (▼) button to adjust height. Press 'Set' to confirm.

Set Age - Use the (▲) or (▼) button to select your age. Press 'Set' to confirm.

Set Gender – Use the (▲) or (▼) button to select female or male icon. Press 'Set' to confirm.

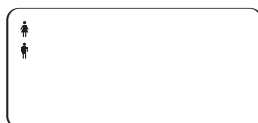Set Fitness Level – Use the (▲) or (▼) button to select fitness level (see Page 4 for definitions). Press 'Set' to confirm.

User memory storage is complete.

User information will show 'Memory' location, height, age and fitness level for 2 seconds each, before unit turns off automatically.

Repeat process for each additional user.

To Clear Memory

To enter new data for a User, complete steps in **Programming Personal Data into Memory** on pages 5 & 6. The new data will override the stored data.


Instructions For Use

Taking Your Body Fat, Body Water, BMR & Muscle Mass Reading:

You must have bare feet to take this measurement. Slightly dampen bottom of feet before stepping on analyzer for best accuracy.

Press the 'User' button to select your 'Memory' location and recall personal information.

User number, height, age and fitness level are each displayed for 2 seconds.

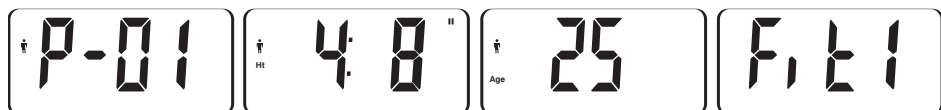(1) Wait for unit to display "0.0"


(2) Step on HealthStation with bare feet. Stand still while your weight is calculated. The weight will be displayed for 3 seconds. Remain standing on the analyzer.


Measurement of body fat, body water, BMR and muscle mass begins.

(3) 'bF%' displays for 2 seconds. Your body fat percentage follows, and is shown for 3 seconds, with the word 'LOW', 'NORMAL', 'MODERATE', or 'HIGH' illuminated. See *Body Fat Percentage Range Chart* on page 2 for recommended body fat levels.

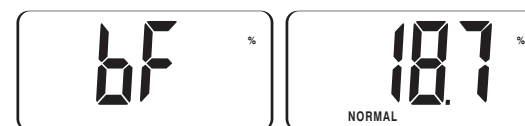(4) 'H2O%' displays for 2 seconds. Your body water percentage follows, and is shown for 3 seconds. See *Body Water Chart* on page 3 for recommended body water levels.

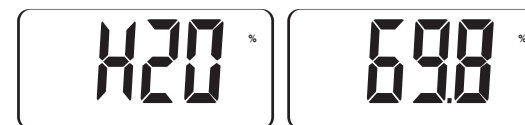(5) 'br' displays for 2 seconds. Your Basal Metabolic Rate follows, and is shown for 3 seconds.

(6) 'M%' displays for 2 seconds. Your muscle mass percentage follows, and is shown for 3 seconds.

Unit turns off automatically.

Troubleshooting

- You must have bare feet when using body fat, total body water, BMR and muscle mass measurements.
- To get the most accurate and consistent reading, wipe your feet with a damp cloth, leaving them slightly moist before stepping on your HealthStation.
- Be sure your feet are placed on the electrodes.

'Err' is displayed

- Body fat is lower than 5% or over 60%
- Improper feet contact when user has socks or shoes on.
- Hard or dry skin can prevent you making proper contact with electrodes on the surface of the analyzer. Try applying a small amount of moisturizer to your feet before using the unit.
- HealthStation is overloaded.


"Lo" is displayed

This is the low battery warning. Battery needs to be replaced.

Battery Replacement

- Open battery cover
- Pull battery out from right side gap in the battery compartment (Picture 1)
- Insert new battery into the left shelter of the battery compartment with "+" facing up (Picture 2)

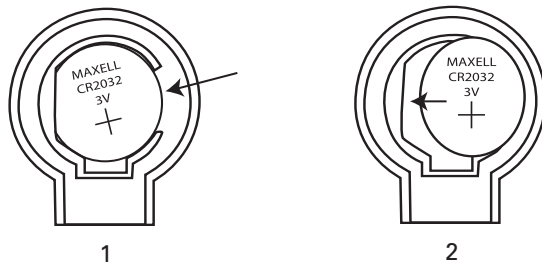- Close battery cover

Care and Maintenance


- Do not disassemble the unit, other than replacing the battery; it contains no user serviceable parts. Damage to the unit due to improper handling will void your warranty.
- Remove battery when the product is not used for a prolonged period of time.
- Clean the unit after use with a lightly dampened cloth. Do not use solvents or immerse the unit in water.
- Avoid excessive impact or vibration, such as dropping it onto the floor.

This product is intended for home/consumer use only; it is not intended for professional use in hospitals or medical facilities.

Note: If you are having any problems operating this HealthStation, please call 800-466-3342, between 8:30am to 5:00pm EST

Specifications

- Uses BIA (Bioelectrical Impedance Analysis) to determine body fat, total body water and muscle mass percentage
- High precision STRAIN GAUGE technology to determine weight
- Strong tempered safety glass platform
- 1.3" Digital display
- Weight Range: 20 to 350 Lb (10 to 160 Kg)
- Weight Graduation: 0.2 Lb (0.1 Kg)
- Body Fat Graduation: 0.1%
- Body Water Graduation 0.1%
- Muscle Mass Graduation: 0.1%
- BMR Graduation: 1 Calorie
- Height Range: 3' 3.4" to 7' 2" (100cm - 220cm)
- Age Range: 10 to 100
- Weight-only feature
- Auto on/off function for convenience and conserving power
- Data can be stored for up to 10 individuals
- User may select unit of measure: pounds (Lb/in), kilograms (Kg/cm), or stones (St/in)
- Low battery indicator
- Operates on one 3 V Lithium battery (included)

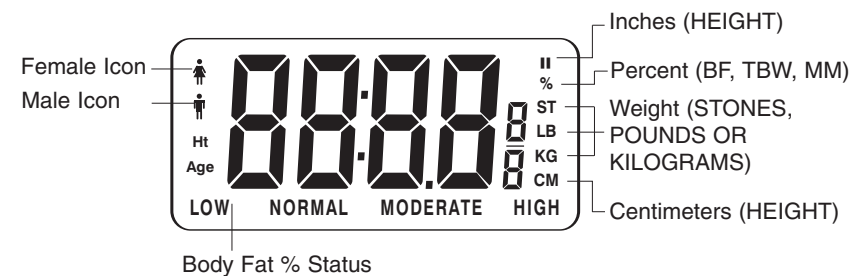LIFETIME WARRANTY

HoMedics, Inc., guarantees this product free from defects in material and workmanship for the life of the original purchaser from the date of original purchase, except as noted below.

This HoMedics product warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; or any other conditions whatsoever that are beyond the control of HoMedics. This warranty is effective only if the product is purchased and operated in the USA. A product that requires modification or adaptation to enable it to operate in any country other than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under warranty. HoMedics shall not be responsible for any type of incidental, consequential or special damages.

For service or repair - do not return this unit to the retailer. Contact HoMedics Consumer Relations at 1-800-HOMEDIC (1-800-466-3342).

Mail To:

HoMedics
Service Center Dept. 168
3000 Pontiac Trail
Commerce Township, MI
48390

e-mail:

cservice@homedics.com

To obtain warranty service on your HoMedics product, either hand deliver or mail the unit and your dated sales receipt (as proof of purchase), postpaid, along with check or money order in the amount of \$5.00 payable to HoMedics, Inc. to cover handling.

Upon receipt, HoMedics will repair or replace, as appropriate, your product and return it to you, postpaid. If it is appropriate to replace your product, HoMedics will replace the product with the same product or a comparable product at HoMedics' option. Warranty is solely through HoMedics Service Center. Service of this product by anyone other than HoMedics Service Center voids warranty.

This warranty provides you with specific legal rights. You may have additional rights which may vary from state to state. Because of individual state regulations, some of the above limitations and exclusions may not apply to you.

For more information regarding our product line in the USA, please visit: www.homedics.com

© 2005 HoMedics, Inc. and its affiliated companies, all rights reserved. HoMedics[®] is a registered trademark of HoMedics, Inc. and its affiliated companies. The beauty of living well.™, TRI-FITNESS™, Know exactly what you're made of.™ and HealthStation™ are trademarks of HoMedics, Inc. and its affiliated companies. All rights reserved.

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